




ONE® Deep Well Full-Cut Punch

Suitable for PK, DALK & DSAEK
for donor graft preparation

Control

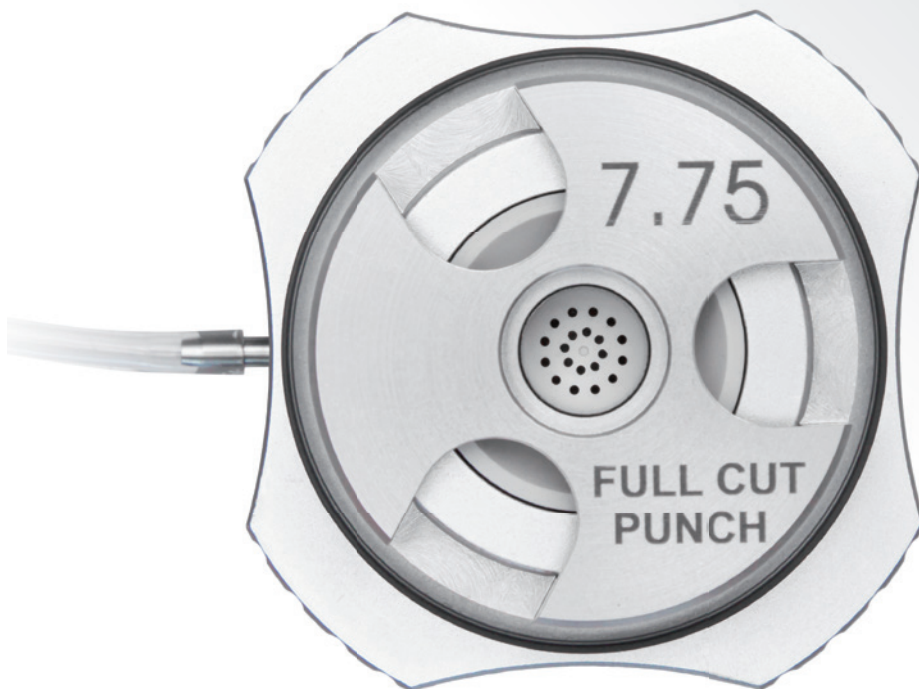

Precision

Ergonomic

ONE STEP FURTHER TO PERFORM YOUR GRAFT PREPARATION TECHNIQUE SUITABLE FOR PK, DALK & DSAEK DONOR GRAFT PREPARATION

The **ONE® Deep Well Full-Cut Punch** is a single-use instrument used for donor cornea. It provides accuracy by using an anatomical well form that allows to center and envelop the cornea.

The syringe creates the vacuum to prevent the cornea from slipping, and ensures that donor cornea stays in place when cutting and removing the blade.

To make sure of a complete 360° full-thickness incision, the blades used are thin and ultra sharp.

In addition, surgeons can control through the three windows.

Full sterile kit, provided with syringe, top blade and base.

**Moria, your best ally to prepare
FULL Clear Uniform Transfixiant graft**

THE DEEP WELL FULL-CUT PUNCH

GOOD CENTERING


Anatomical deep well with double curvature to envelop the cornea in its most natural position

PRECISE

20 holes distributing vacuum for suction and secure graft

ERGONOMIC

Comfortable & stable

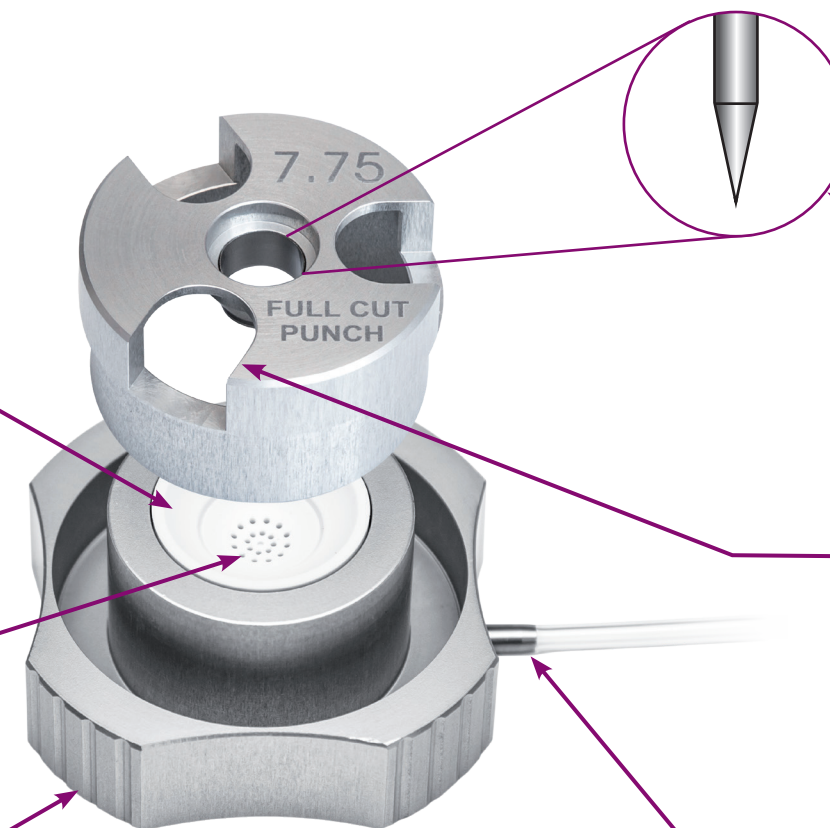

DOUBLE-BEVELLED SYMETRICAL BLADE

for a vertical uniform incision

GOOD VISIBILITY & GENTLE CONTROL

through 3 lateral windows


SECURE

Good suction to hold donor cornea in place

REFERENCES

14 sizes from 6.00 mm to 10.00 mm

Deep Well Full-Cut Punch 17214DXXX	
References	Incision diameters
17214D600	6.00 mm
17214D650	6.50 mm
17214D675	6.75 mm
17214D700	7.00 mm
17214D725	7.25 mm
17214D750	7.50 mm
17214D775	7.75 mm
17214D800	8.00 mm
17214D825	8.25 mm
17214D850	8.50 mm
17214D900	9.00 mm
17214D950	9.50 mm
17214D1000	10.00 mm

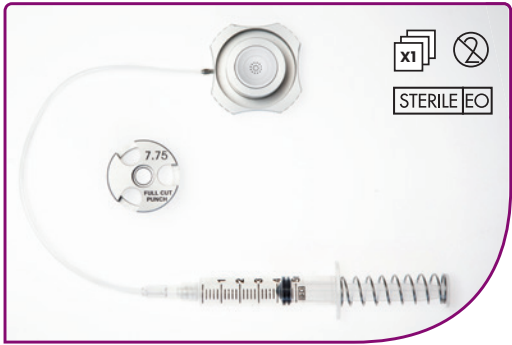

FDA registered & CE marked



MORIA #67086EN-B-07.2024

MORIA SA
27, rue du Pied de Fourche
03160 Bourbon L'Archambault
FRANCE
Phone: +33 (0)1 46 74 46 74
Fax: +33 (0)1 46 74 46 00
moria@moria-int.com
www.moria-surgical.com

MORIA Inc
1050 Cross Keys Drive
Doylestown, PA 18902
USA
Phone: (800) 441 1314
Fax: +1 (215) 230 7670
orders@moriausa.com
www.moria-surgical.com

Moria

Moria Japan K.K.
Arcadia Building 6F
1-12-3 Kanda
SudachoChiyoda-Ku
Tokyo 101-0041
JAPAN
Phone: 81-3-6260-8309
Fax: 81-3-6260-8310
moria@moriajapan.com
www.moriajapan.com

Moria COMMERCIAL (CHINA) CO., LTD.
上海目利亚贸易有限公司
Room H 6F Kaili Building NO.432
West Huai Hai Road, Changning district, Shanghai, 200052, P.R.C
中国上海市长宁区淮海西路432号
凯利大厦6楼08室
Phone/Fax : +86 021 52586095
moriachina@moria-int.com