



# Keratoplasty

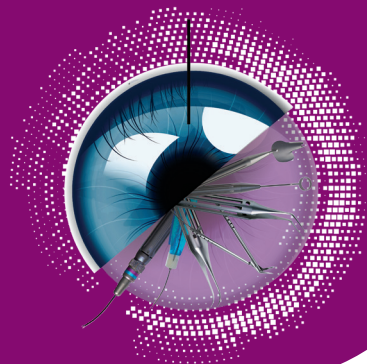

**ONE® Deep Well Full-Cut Punch**  
Suitable for PK, DALK & DSAEK  
donor graft preparation

*One step further to perform your graft preparation*

**Control**

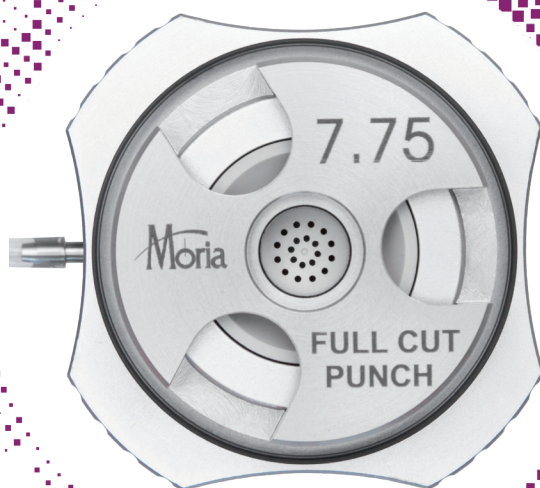

**Precision**

**Ergonomic**





# One step further to perform your graft preparation technique. Suitable for PK, DALK & DSAEK donor graft preparation.



The **One® Deep Well Full-Cut Punch** is a single-use instrument used for donor cornea. It provides accuracy by using an anatomical well form that allows to center and envelop the cornea.

The syringe creates the vacuum to prevent the cornea from slipping and ensures that cornea stays in place when cutting and removing the blade.

To avoid undercutting, the blades used are thin and ultra sharp.

In addition, surgeons can control through the three windows.

Full sterile kit, provided with syringe, top blade and base.

**Moria, your best ally to have**  
**FULL Clear Uniform Transfixiant graft preparation**

## One® Deep Well Full-Cut Punch

### ■ Good centering

Anatomical deep well with double curvature to envelop the cornea in its most natural position

### ■ Precise

20 holes distributing vacuum for suction and secure graft

### ■ ERGONOMIC

Comfortable & stable

### ■ Double-bevelled symetrical blade ■

for a vertical uniform cut

### ■ Good visibility & gentle control ■

through 3 lateral windows

### ■ Secure ■


Good suction to hold donor cornea in place



# References

14 sizes from 6.00 mm to 10.00 mm  
with an increment of 0.25 mm between diameter 6.50 mm and 9.00 mm

One® Deep Well Full-Cut Punch	
References	Sizing diameter
17214D600	6.00 mm
17214D650	6.50 mm
17214D675	6.75 mm
17214D700	7.00 mm
17214D725	7.25 mm
17214D750	7.50 mm
17214D775	7.75 mm
17214D800	8.00 mm
17214D825	8.25 mm
17214D850	8.50 mm
17214D875	8.75 mm
17214D900	9.00 mm
17214D950	9.50 mm
17214D1000	10.00 mm

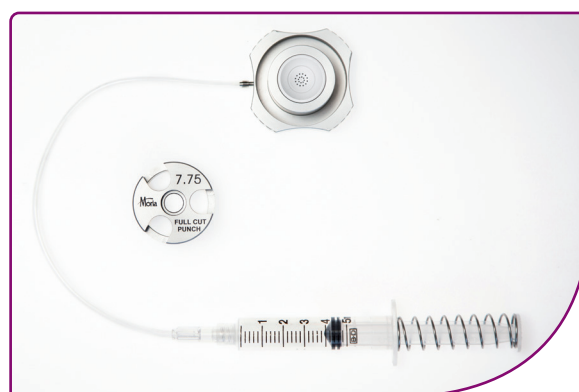

FDA registered




To obtain more information

**MORIA SA**  
15, rue Georges Besse  
92160 Antony  
FRANCE  
Phone: +33 (0) 1 46 74 46 74  
Fax: +33 (0) 1 46 74 46 00  
[www.moria-surgical.com](http://www.moria-surgical.com)



**MORIA Inc**  
1050 Cross Keys Drive  
Doylestown, PA 18902  
USA  
Phone: (800) 441 1314  
Fax: +1 (215) 230 7670  
[www.moria-surgical.com](http://www.moria-surgical.com)

MADE IN FRANCE

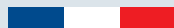

MORIA #67086EN-A-02.2022

**Moria**  
Ophthalmic Instruments